



Universal Connector Range

Products available with tubing featuring a universal connector



Oxygen and Aerosol Therapy


Universal Connector

The products included in this information sheet are available with tubing featuring a universal connector, which can be attached to a 6 mm flow meter spigot or fitted directly onto a 9/16-18 UNF Diameter Index Safety System (D.I.S.S.) flow meter.



Medium concentration oxygen mask

The Intersurgical EcoLite™ medium concentration oxygen mask is lightweight, and comfortable. Every aspect of the design has been addressed, improving levels of comfort and fit. Not only comfortable for the patient but also comfortable for the environment. The elimination of PVC from the Intersurgical EcoLite masks has lowered their environmental impact.

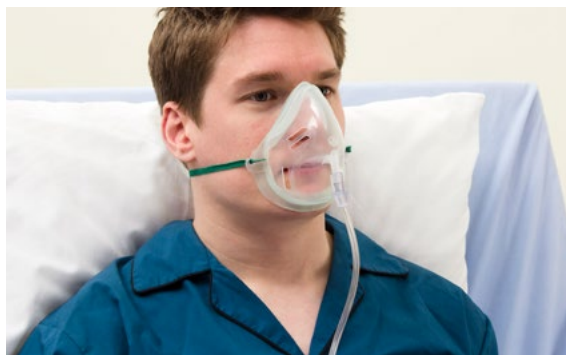

Code	Description	Tube length	Box Qty.	
1135017	Intersurgical EcoLite, adult, medium concentration oxygen mask and tube with universal connector	2.1 m	30	
1196019	Intersurgical EcoLite, paediatric, medium concentration oxygen mask and tube with universal connector	2.1 m	40	

[Make an enquiry](#)

[Watch the video](#)

High concentration oxygen mask

The Intersurgical EcoLite adult high concentration mask is designed to be comfortable for the patient, with a soft face seal. The product features two integrated safety valves, which allow the patient to breathe room air in the event oxygen supply failure.



Code	Description	Tube length	Box Qty.	
1181017	Intersurgical EcoLite, adult, high concentration oxygen mask and tube with universal connector	2.1 m	24	

[Make an enquiry](#)

[Watch the video](#)



Silente™ venturi mask kit

The Silente has been designed to operate with minimal noise output, with the large diameter body providing space for clear, widely spaced calibration and incorporating a “pull-twist” action to help minimise the risk of accidental adjustment.

The Silente venturi valve has been designed to deliver fixed concentrations of oxygen, that may be adjusted to suit clinical requirements.



Code	Description	Tube length	Box Qty.	
1137016	Silente, Intersurgical EcoLite™, adult, adjustable venturi mask kit and tube with universal connector	1.8 m	30	

[Make an enquiry](#)
[Watch the video](#)

Sentri™ ETCO₂ mask

The Sentri Mask is an adult, medium concentration Intersurgical EcoLite™ oxygen mask with a ETCO₂ monitoring line connector enabling a secure luer lock connection, the other end of which is connected to a capnometer. The Sentri mask is used for measuring end-tidal CO₂ (ETCO₂) in expired gases while supplying oxygen to the patient.



Code	Description	Tube length	Box Qty.	
1143018	Sentri Intersurgical EcoLite, adult, mask with CO ₂ monitoring line, filter and tube with universal connector	2.1 m	30	

[Make an enquiry](#)

Cirrus™2 nebuliser mask kit

The Cirrus2 nebuliser is designed to deliver medication for tracheobronchial deposition. Nebuliser performance specifications are detailed in documentation provided with the products, and are available upon request.

The Cirrus2 mask kits incorporate the Intersurgical EcoLite aerosol mask, this non-PVC mask reduces the environmental impact of this product.



Code	Description	Tube length	Box Qty.	
1453014	Cirrus2 nebuliser, adult, Intersurgical EcoLite mask kit and tube with universal connector	2.1 m	36	
1454017	Cirrus2 nebuliser, paediatric, Intersurgical EcoLite mask kit and tube with universal connector	2.1 m	40	
1455002	Cirrus2 nebuliser, mouthpiece kit and tube with universal connector	2.1 m	35	

[Make an enquiry](#)
[Watch the video](#)

HOT Top™2 nebulisers

The HOT Top2 nebuliser is an open nebuliser incorporating twin entrainment ports for high output flows resulting in reduced treatment times and improved patient compliance.

In addition, the HOT Top2 nebuliser includes a graduated nebuliser cup which provides guidance on how much medication is left in the nebuliser and a quick-release fitting which enables the nebuliser to be attached and detached with a quarter-turn. This speeds up the filling process and provides a reliable seal between the cup and the nebuliser body.



Code	Description	Tube length	Box Qty.
1463016	HOT Top2 nebuliser, adult, Intersurgical EcoLite mask kit and tube with universal connector	1.8 m	30
1468016	HOT Top2 nebuliser, paediatric, Intersurgical EcoLite mask kit and tube with universal connector	1.8 m	36

[Make an enquiry](#)
[Watch the video](#)

Oxygen tube with universal connector

Oxygen tubing featuring a universal connector, which can be attached to a 6 mm flow meter spigot or fitted directly onto a 9/16-18 UNF Diameter Index Safety System (D.I.S.S.) flow meter in a 2.1m length, the tube is kink resistant helping to prevent occlusion.



Code	Description	Tube length	Box Qty.
1174021	Oxygen tube with universal connector	2.1 m	100

Note: Not all products are available in all markets, please contact your local office for details.

IS10.41 • Issue 2 03.26



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

UK • Australia • Benelux • Brazil • Canada • China • Colombia • Czech Republic • Denmark • France • Germany • India • Ireland • Italy • Japan • Lithuania • Malaysia • Mexico • Norway • Philippines • Portugal • Russia • South Africa • Spain • Sweden • Taiwan • Turkey • USA