




STARMed CaStar

hood for CPAP therapy without manometer

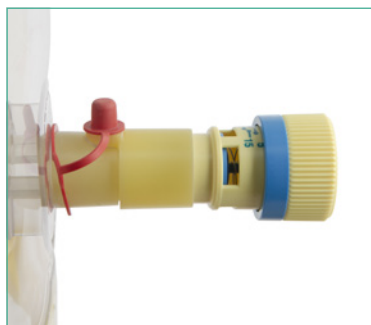

Critical Care ▪ Patient Interfaces ▪ **STARMed**




STARMed CaStar hood for CPAP therapy without manometer

The CaStar is a comfortable, versatile and lightweight CPAP patient interface designed to offer an alternative to traditional methods of CPAP delivery. The CaStar is well tolerated for extended periods and can be used in semi-recumbent, supine and prone positions. StarMed CPAP hoods feature a built-in manometer to check the presence of pressure within the hood. The manometer working principle is based on a minimal intentional leak. To further reduce the occurrence of aerosol generating procedures in the clinical environment, this feature has been removed from this CPAP hood and an external manometer accessory (CP929) can be applied.

Technical Features

- Manufactured from phthalate-free materials
- 22M/22F connectors compliant to ISO standard
- Comfortable and atraumatic underarm fastening straps
- Two versatile sealed access ports suitable for single or multiple lumen probes or catheters
- Patient access port, fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

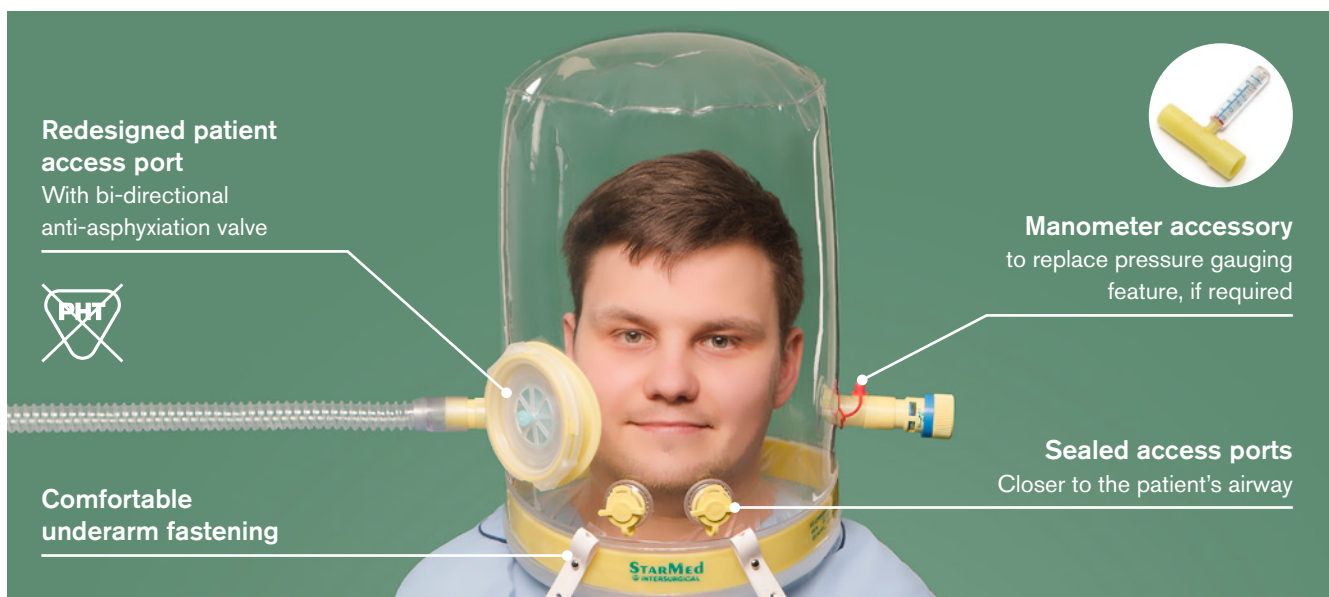

CaStar CPAP hood with bi-directional anti-asphyxiation valve and patient access port, without manometer*

Code	Size	Neck circumference	Box Qty.
CP041XXL/2	Extra Extra Large	≥ 50 cm	5
CP041XL/2	Extra Large	45–52 cm	5
CP041L/2	Large	40–47 cm	5
CP041M/2	Medium	34–41 cm	5
CP041S/2	Small	27–34 cm	5
CP041XS/2	Extra Small	17–27 cm	5
CP929	22M/22F adaptor with built-in manometer		5

[Make an enquiry](#)
[Find out more](#)
[▶ Watch the video](#)

*PEEP valve, tubing and manometer accessory not included and have to be ordered separately. Please refer to the website for the range of accessories.



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



Manufacturer:
Intersurgical S.p.A.
Mirandola (MO) Italy

Please think before you print

Save energy and paper. If you must print this information sheet please print it double sided.

IS1.68_INT • Issue 1 04.25