




One-piece Guedel airways

Safer by design

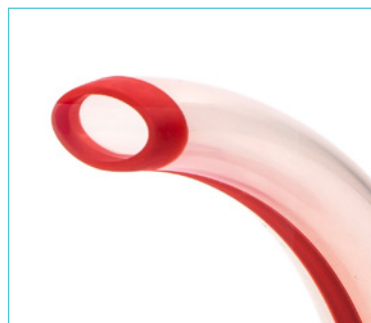

One-piece Guedel airways

Safer by design

By using the latest manufacturing technology Intersurgical has produced the first major breakthrough in Guedel airway design for thirty years.

Combining two materials with very different physical characteristics into one moulding, our One-piece design eliminates the separate bite block found in all conventional Guedel airways.

Loose or detached bite blocks constitute a real threat to patient safety and have been the subject of several hazard notices in the past. Our One-piece design overcomes all of the associated risks.

As well as improved safety, the One-piece Guedel airway also incorporates a soft tip reducing the risk of trauma to the patient. We also offer a sterile option, providing a further choice.

Rigid body
To maintain a patent airway

One-piece design
With no separate inner bite block – safer by design

Nine colour-coded sizes
To meet all patient requirements

Soft tip
Reduces risk of tissue damage

Softer second material
Reduces the potential for dental damage

Clear colour-coding
For easy identification





ISO 5364:2016 – Colour coding guidance

In the most recent version of ISO 5364:2016 Anaesthetic and Respiratory Equipment – Oropharyngeal Airways, a new colour coding requirement was implemented, aiming to improve patient safety and provide consistency across products from different manufacturers.

In order to comply with the colour coding system specified in ISO 5364:2016, it was necessary for Intersurgical to update four sizes of Guedels. The sizes that have changed are as follows:

- ISO size 9.0; from Orange to Yellow
- ISO size 7.0; from Yellow to White
- ISO size 6.5; from White to Brown
- ISO size 3.5; from Pink to Light Green



The four sizes of Intersurgical Guedel airways that have had their colours changed to conform to the updated standard.

Each department will need to consider locally how this change will be implemented as it may take time for existing stock to be used.

Product code	ISO size
1115120 (S*)	12.0
1114100 (S*)	10.0
1113090 (S*)	9.0
1112080 (S*)	8.0
1111570 (S*)	7.0
1111065 (S*)	6.5
1110055 (S*)	5.5
1100050 (S*)	5.0
1000035 (S*)	3.5

Previous Intersurgical colour	
Purple	
Red	
Orange	
Green	
Yellow	
White	
Grey	
Blue	
Pink	

New ISO compliant colour	
Purple	
Red	
Yellow	NEW
Green	
White	NEW
Brown	NEW
Grey	
Blue	
Light Green	NEW

This colour change effects all of Intersurgical's Guedel Airway products, including bulk-packed options.

Sterile option available

(S*) Add an S to the seven digit code number for the sterile version of this product eg. 5000000S

Softer tip

By continuing to prioritise patient safety, Intersurgical has further increased the softness of the tip of all Guedel airways.

Guedel airways that utilise hard and rigid materials at the tip represent a significant threat of causing trauma when inserted. By using an even softer material to form the tip of the One-piece Guedel airway, the potential to cause trauma to the patient's airway upon insertion is further reduced.



Code	Description	Colour	Size	ISO	Box Qty.
1115120 (S*)	One-piece Guedel airway		5	12.0	70 (70*)
1114100 (S*)	One-piece Guedel airway		4	10.0	90 (80*)
1113090 (S*)	One-piece Guedel airway		3	9.0	90 (80*)
1112080 (S*)	One-piece Guedel airway		2	8.0	100 (100*)
1111570 (S*)	One-piece Guedel airway		1.5	7.0	100 (100*)
1111065 (S*)	One-piece Guedel airway		1	6.5	20 (20*)
1110055 (S*)	One-piece Guedel airway		0	5.5	25 (25*)
1100050 (S*)	One-piece Guedel airway		00	5.0	25 (25*)
1000035 (S*)	One-piece Guedel airway		000	3.5	25 (25*)

- [Make an enquiry](#)
- [Find out more](#)
- [Watch the video](#)

Lower environmental impact product Sterile option available

(S*) Add an S to the seven digit code number for the sterile version of this product eg. 5000000S

IS6.10 ▪ Issue 7 07.24



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

Please think before you print
 Save energy and paper.
 If you must print this information sheet please print it double sided.